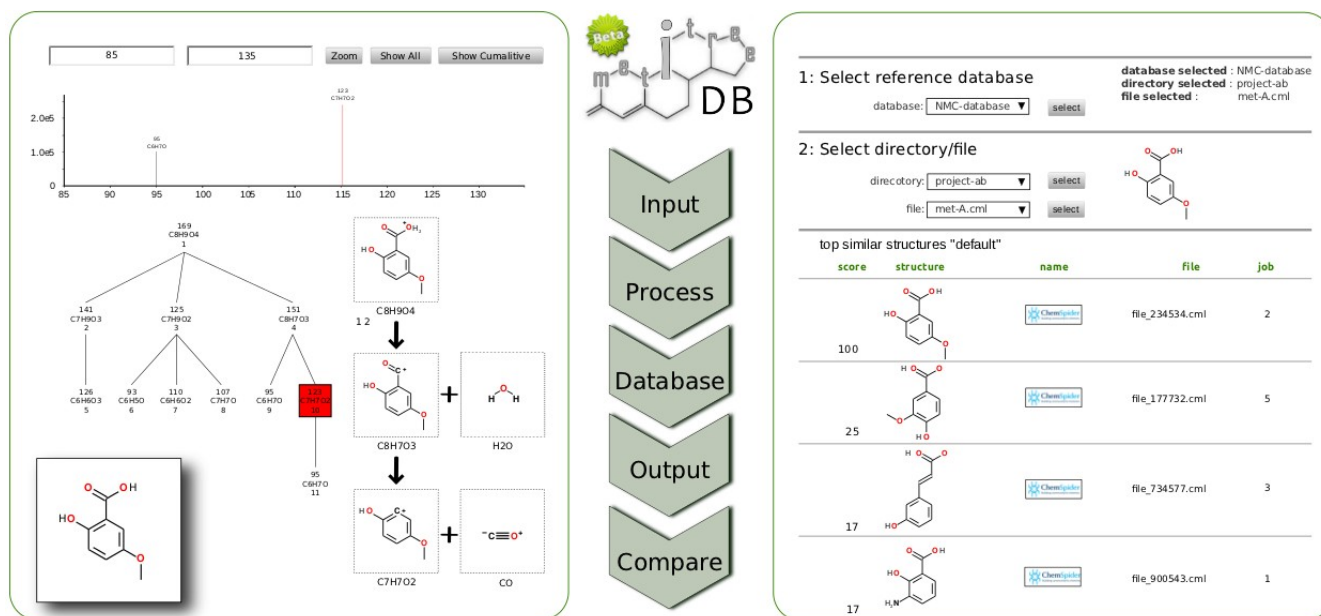


# METITREE USER GUIDE

version 0.1  
First edition: June 2012

This User Guide provides an overview of several functional features of the MetiTree web application. Use the content list to enter detailed information on how to use MetiTree.



MetiTree is a user-friendly, web application dedicated to organize, process, share, visualize, and compare high resolution multi-stage mass spectrometry ( $MS^n$ ) data. It integrates dedicated features to export and visualize complex  $MS^n$  data, facilitating the exploration and interpretation of metabolomics experiments.

## Table of Contents

Upload $MS^n$ data: .....	2
Process and visualize .....	3
Storage and manage database(s) with $MS^n$ data .....	6
Querying an existing database(s) for similar $MS^n$ data: .....	8
Visualize $MS^n$ data with chemical structure of the fragments: .....	9
FAQs: .....	10
Biography: .....	10

## Upload MS<sup>n</sup> data

Upload files with MS<sup>n</sup> mass spectrometry data and store them into data directories.


MetiTree accepts right now the uploading of mzXML files (version 3.1) containing MS<sup>n</sup> mass spectrometry data. MetiTree allows the creation of directories for grouping mzXML files, assisting the organization of the data according projects or topics.

1. Select “**MS<sup>n</sup> data**” (tab), to open the “**MS<sup>n</sup> data**” page.
2. From the '**MS<sup>n</sup> data**' you can perform the functions listed below:
  - Create a new directory
  - Modify an existing directory
3. To create a new directory, first provide a name for the new directory. Click '**create**' to save your directory.




### create a directory

Use directories to group files. Directories can also be used to process multiple files at once.

name :

4. To modify an existing directory, select one of the directories  by clicking the magnifying glass.

### directories

	name	file	created
	mi-files	1	2012-01-16
	publication-files	20	2012-01-16
	testFilesMichael	1	2012-01-16

5. The page “**Files**” will be displayed. The files added to the directory will be listed.

files in: *metabolites* (change directory / add files)

				sort by: name / date							
ID	file	name		pubchem	chemspider	owner	delete				
 19	dilu_c.mzxml	Glutathione				you					
 18	dilu_d.mzxml	Glutathione				you					
 17	hmdb_a.mzxml	6-hydroxy-m-Anisic acid				you					
 16	lab_c.mzxml	5-Hydroxylysine				you					
 15	hmdb_f.mzxml	N-Acetylserotonin				you					
 14	hmdb_g.mzxml	3-(2-Hydroxyethyl)-1H-indol-5-ol				you					

6. To see extended information of a file, click magnifying glass. 

files in: *metabolites* (change directory / add files)

sort by: name / date					pubchem	chemspider	owner	delete
ID	file	name						
19	dilu_c.mzxml	Glutathione					you	
18	dilu_d.mzxml	Glutathione					you	

2d image	InChI	comments
	<div>InChI=1S/C10H17N3O6S/c11-5(10(18)19)1-2-7(14)13-6(4-20)9(17)12-3-8(15)16/h5-6,20H,1-4,11H2,(H,12,17)(H,13,14)(H,15,16)(H,18,19)/t5-,6-/m0/s1</div> <div>update/save</div>	<div></div> <div>update/save</div>

7. Select “**add files**” to add new files from your desktop to this directory. Finally, select “**submit**” to save it.

You can upload a single mzXML file or a zip-file with multiple mzXML files. mzXML file is an open data format for storage and exchange of mass spectroscopy data, developed at the SPC/Institute for Systems Biology. All known converters to convert vendor specific files to mzXML files are summarized in the page (<http://tools.proteomecenter.org/wiki/index.php?title=Formats:mzXML>)

8. Once the file is uploaded you eventually can specify the InChI identifier of the neutral compound. To save the InChI identifier you have to select the “**update**” button. Automatically, the file will be cross-linked to the pubchem and chemspider databases.
9. The files can be sorted by name or data created by selecting “**name/date**”.

## Process and visualize

MS<sup>n</sup> data:

MetiTree integrates the MEF tool (Rojas-Chertó et al., 2011) and uses its functionality to process MS<sup>n</sup> data. The required input to process MS<sup>n</sup> data are mzXML files (<http://en.wikipedia.org/wiki/MzXML>) and the settings of the processing parameters.

Processing parameters are grouped into those to extract the mass spectrometry information (m/z, intensity, and scan number) and those to enrich the MS data with chemical information (chemical elements and number of atoms). MetiTree allows the processing of individual or multiple files at the same time. Furthermore, the same mzXML file can be processed several times with different sets of processing parameters. Results and parameters information are stored to allow posterior revision.

Once the data is processed, it can be displayed using the spectral tree viewer. Here, the user

can select a peak or a node, which represents a fragment. For each selection, the corresponding mass spectrum is displayed, together with the reactions that connect the parent ion with the selected fragment. The structures of the fragments are only visible when they have been assigned. The results generated by MetiTree can be exported to different file formats for further analysis. The CML format allows managing complex chemical content that other software can import to process it even further. The PDF format can be used to visualize the MS<sup>n</sup> data, which makes it suitable to export images into reports or publications.

1. On the navigation bar select “**Process**” (tab)
2. To process a file you must first select the directory where the file is found. To process all files of a determined directory select “**process directory**”.

select what to process :

---

**please select a directory to process or select files from**

publication-files ▼

select file(s) from directory to process or process directory

3. If you want to process specific files from the directory, select “**select file(s) from directory to process**”. A list of files will be displayed. Check those files to process and select “**process files**”. If you like to select all files at once press the box with the sign plus.

**files in directory *publication-files***

Process files

	file	created
<input checked="" type="checkbox"/>	HMDB_a.mzXML	2012-01-15
<input checked="" type="checkbox"/>	HMDB_b.mzXML	2012-01-15
<input checked="" type="checkbox"/>	HMDB_c.mzXML	2012-01-15
<input type="checkbox"/>	HMDB_d.mzXML	2012-01-15
<input type="checkbox"/>	HMDB_e.mzXML	2012-01-15

4. The “**Define settings**” page will be displayed introducing the following processing settings:
  - MZ gap (bin size) = Minimal distance between adjacent peaks. Peaks with smaller intensity being excluded.
  - Signal to noise threshold = Signal to noise threshold allowed during the peak picking.
  - Elements = It defines upper and lower number of atoms admitted for each element in the elemental composition of the ion (e.g. C1..15, H1..9, O0..4).
  - Rules = It defines the constraint-rules that are applied to verify the elemental compositions generated. It distinguishes between the Nitrogen Rules and the RDB rules.
  - Mass accuracy settings = Mass accuracy error allowed per MS level.

Please change the settings or continue with the defaults

### Peak picking settings

MZ gap (bin size) :  (5 e.g 0.5 m/z)

Signal to noise threshold :

### Filter settings

Elements :  ('C1..50,H1..100,N0..10,O0..10' or 'C10,N4,H2')

Rules :

☒ nitrogenR

☒ RDBER

### Mass accuracy settings

MS Level 1 :

MS Level 2 :

MS Level 3 :

MS Level 4 :










MS Level 5 :

**Process**

5. To submit the process select “**Process**”.

Status JOB#000000000022 (...finished!)

#### Output

name	compound	mzxml 	cml 	export (csv/tab) 	visuals 	created
dilu_d.mzXML	Glutathione 					2012-06-15

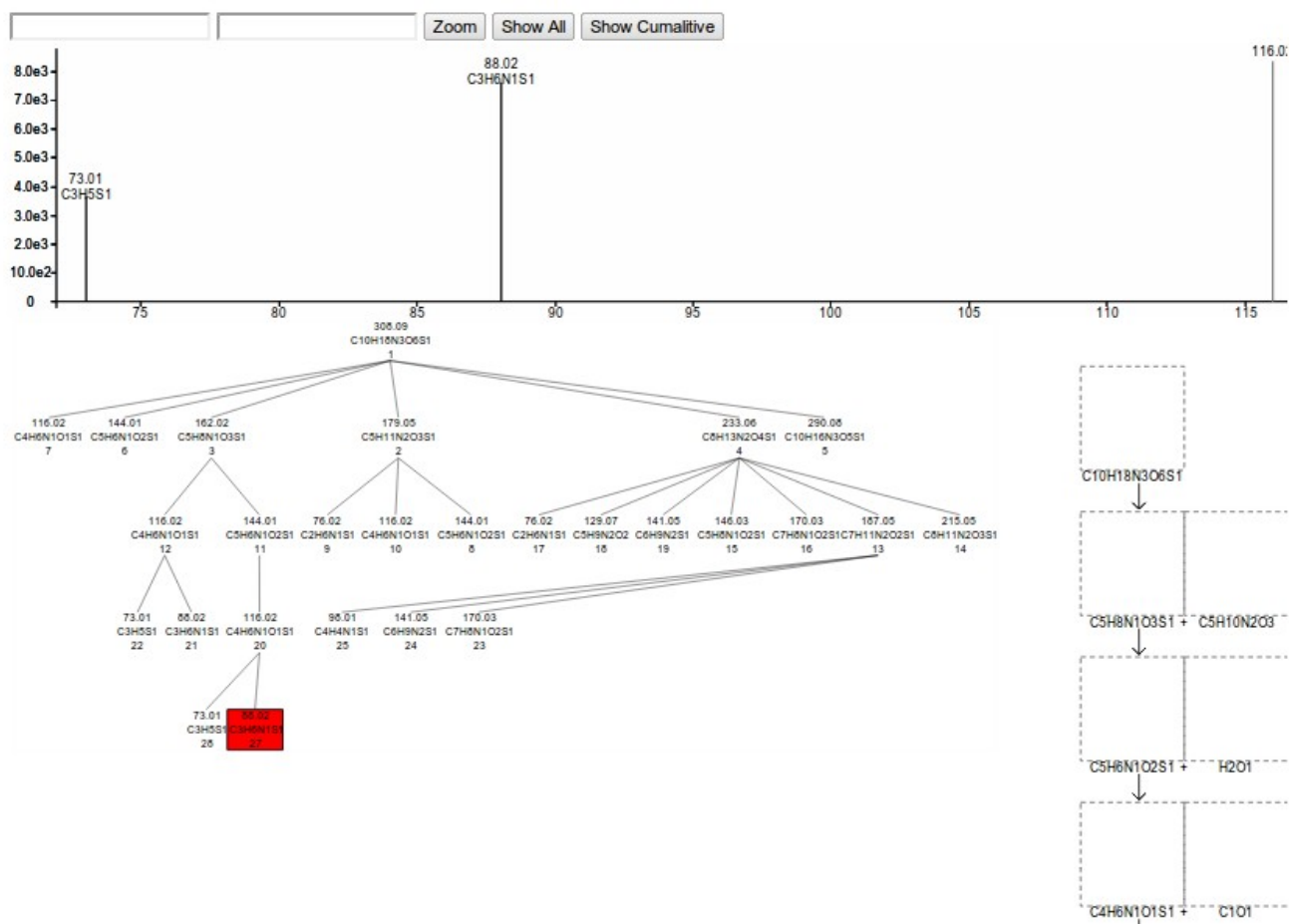
#### Settings

setting	value
accuracy-L1	15
accuracy-L2	15
accuracy-L3	15
accuracy-L4	15
accuracy-L5	15
elements	O6..6,N3..3,H18..18,C10..10,S1..1
mzgap	0.5
rulenitrogenr	nitrogenR
rulerdber	RDBER
snthresh	1.0

#### Job Log

Fri Jun 15 17:19:46 CEST 2012 - Job queued  
Fri Jun 15 17:19:46 CEST 2012 - Job started  
Fri Jun 15 17:20:32 CEST 2012 - **ERROR:** Creation of PDF formatted files failed on file: dilu\_d.mzXML.pdf...

6. The “**Process**” page will be displayed containing information about the progress of the different processing jobs.
7. It is possible to visualize the different results.
  - mzXML = original uploaded mzXML file used to process the spectral data.
  - cml = chemical markup language file enriched with chemical information.
  - trees = Tree Viewer interconnects three MS<sup>n</sup> items: the spectrum, which contains mass peaks, the fragmentation tree, which contains fragment nodes/elemental compositions, and the fragmentation reactions, which contain structures.
  - csv = comma separated values file with enriched chemical information
  - pdf = portable document format to visualize the spectral data.
8. Select “**master tree**” to visualize spectral data with the tree viewer



## Storage and manage database(s) with MS<sup>n</sup> data


Processed MS<sup>n</sup> data can be stored in one or multiple internal databases. Because the users




are organized in groups, they can share files and libraries with other group members. Only the administrator user has the faculty to create a new group and assign users to it.

1. On the navigation bar select the “**Databases**” (tab)
2. The '**Databases**' page will be displayed. From this page, you have access to the functions listed below:
  - Create a new database
  - Modify an existing database
3. To create new database, provide a name for the new database and click '**add**' to save your database.

**add a new database**

name :

4. To modify an existing database, select one of  the magnifying glasses in front of the databases.

databases			
	name	file	created
	publication-database	13	2012-01-15 04:30:46.08
	ronnie-albert	1	2012-01-15 04:31:42.24
	testDatabaseMichael	1	2012-01-15 04:31:56.617

5. Now modify the files to be contained in the database as required.



#### publication-database

filter on job

All jobs

filter

Add files to database





Files that can be added to this database.



























 add all

Files in database

The data in this database originates from the files below.

 remove all

File	Job	Created
 hmdb00434neg_01.mzXML.cml	2	 2012-01-15 04:31:45.663
 hmdb00434neg_01.mzXML.cml	3	 2012-01-15 04:32:01.519

File	Job	Created
 HMDB_a.mzXML.cml	19	 2012-01-15 04:41:43.521
 HMDB_b.mzXML.cml	13	 2012-01-15 04:39:27.985
 HMDB_c.mzXML.cml	15	 2012-01-15 04:40:35.075
 HMDB_d.mzXML.cml	5	 2012-01-15 04:33:15.941
 HMDB_e.mzXML.cml	8	 2012-01-15 04:35:07.403
 HMDB_f.mzXML.cml	17	 2012-01-15 04:41:06.914
 HMDB_g.mzXML.cml	16	 2012-01-15 04:40:55.965
 dilu_d.mzXML.cml	21	 2012-01-15 04:42:15.268
 dilu_e.mzXML.cml	20	 2012-01-15 04:42:11.831
 lab_a.mzXML.cml	11	 2012-01-15 04:37:35.395
 lab_b.mzXML.cml	6	 2012-01-15 04:33:22.836
 lab_c.mzXML.cml	18	 2012-01-15 04:41:31.434
 lab_d.mzXML.cml	14	 2012-01-15 04:39:58.725

6. On the right side you have the files that actually are in the database. On the left side you have those files that you can upload to the database. Through “filter on job” you can filter the files by specific job number.

## Querying an existing database(s) for similar MS<sup>n</sup> data

MetiTree integrates the functionality to query for similar MS<sup>n</sup> data stored in your own personal library. The results are presented in a list showing the structures of the compounds with most similar MS<sup>n</sup> data together with the corresponding similarity value. The results are ranged between 0-100. A value near to 100 indicates that the MS<sup>n</sup> data is highly similar which would be a complete identification; while a value close to 0 illustrates that it is very different (Rojaschertó et al., 2012). We expect MetiTree to contribute two fold in the identification of metabolites. On one hand, we are interested to know if a fragmentation tree of the same compound is present in the library, which would be a partial identification. On the other hand, if similar fragmentation data is found, we can use this knowledge to give hints (e.g. substructure information) which structure the unknown compound could have.

1. Select the “search/compare” (tab), to open the “database search” page.
2. Choose a database to query “database”, and the spectral tree file to be compared in “directory” > “file”. To see the database search result push “search/refresh” button.



1: select reference database:

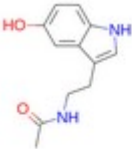
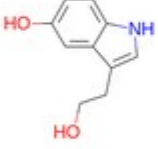
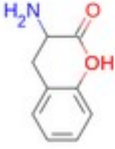
database :

2: select query directory/file:

directory :

file :

3. The '**database search**' page will display the results into a list of compounds ranked from more similar to a less similar. On the left side it will show the similarity value when the fragmentation trees are compared to each other.

score	structure	name	file	job
100.0		PubChem	HMDB_f.mzXML.cml	15
17.51		PubChem	HMDB_g.mzXML.cml	14
6.66		PubChem	HMDB_e.mzXML.cml	6

## Visualize MS<sup>n</sup> data with chemical structure of the fragments

The spectral tree viewer is capable to show the chemical structure of the fragments in the MS<sup>n</sup> data. At the moment it is only possible if you supply it with a JSON file describing the MS<sup>n</sup> data.

1. On the bottom of any page, select the "**MsnViewer Demo**".

2. You can upload an example by selecting **“Load example”** > **“view”**.  
You can modify the json message and add your own values.

Paste a JSON formatted Tree Structure to view with the MSnViewer.js application (load example)

```
{
  "1" : {
    "id" : "1",
    "parent" : "",
    "mass" : "164.01",
    "intensity" : "2601726.5555",
    "inchi" : "InChI=1S/C5H9N03S/c1-3(7)6-4(2-10)5(8)9/h4,10H,2H2,1H3,(H,6,7)
(H,8,9)/p+1,InChI=1S/C3H7N02S/c4-2(1-7)3(5)6/h2,7H,1,4H2,(H,5,6)/p+1",
    "inchiloss" : ""
  },
  "2" : {
    "id" : "2",
    "parent" : "1",
    "mass" : "122.01",
    "intensity" : "378006.0052",
    "inchi" : "InChI=1S/C3H7N02S/c4-2(1-7)3(5)6/h2,7H,1,4H2,(H,5,6)/p+1",
    "inchiloss" : "InChI=1S/C2H20/c1-2-3/h1H2"
  },
  "3" : {
    "id" : "3",
    "parent" : "1",
  }
}
```

View

3. After the selected view run, the MS<sup>n</sup> viewer screen will be shown.

## FAQs

Where can I address for more help?

For any questions that could not be answered through this guide, you can send us an email directly to [info@metitree.nl](mailto:info@metitree.nl). We will answer you as soon as possible.

How many files might I upload to process?

At the moment there is not a limit during the uploading process. If the number is big may be you can collapse the server.

## Biography:

Rojas-Chertó, M., Kasper, P. T., Willighagen, E. L., Vreeken, R., Hankemeier, T., & Reijmers, T. (2011). Elemental Composition determination based on MSn. *Bioinformatics*, 27(17), 2376–2383. Oxford Univ Press. doi:10.1093/bioinformatics/btr409

Rojas-chertó, M., Peironcely, J. E., Kasper, P. T., Van, J. J. J., Hoof, D., Vos, R. C. H. D., Vreeken, R., Hankemeier, T., & Reijmers, T. (2012). Metabolite identification using automated comparison of high resolution MSn spectral trees. *Anal. Chem.* doi:10.1021/ac2034216